Gebrauchsanweisung Art.-Nr. 101763

schwa-medico
MEDIZINTECHNIK

SALIMA Bodywave® KLANGLIEGE
SALIMA Bodywave® KLANGMATTE

Inhalt

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SALIMA Bodywave® SOUND RECLINER
SALIMA Bodywave® SOUND MAT

Content

What You Should Know About The SALIMA Bodywave® Sound Recliner/Mat	15
General	16
Purpose Of Use	16
Safety Instructions	16
Areas Of Application (Indications)	16
Contra-Indications	16
Side Effects	16
Description Of The SALIMA Bodywave®-Sound Recliner/Mat	17
Start-Up	18
Control Unit	18
Power Connection	18
Master Display	18
Program Selection	18
Recommendation	19
Buttons and Controls	19
Application Duration	20
Programs Overview	20
Program 1 - Relaxation	20
Program 2 - Revitalisation	20
Program 3 - Relief For The Legs	20
Program 4 - Lower Abdomen, Lumbar Vertebral Column	20
Program 5 - Neck/Breast	20
Program 6 - Happy Day	20
General Information	21
Care And Cleaning	21
Delivery Content	21
Resettings, Modifications and Repairs	22
Warranty	22
Guarantee	22
Technical Data	22
Explanation Of The Illustrations	23

Dear Customer,

We are very pleased you have decided to purchase a SALIMA Bodywave® Sound Recliner/Mat. These operating instructions will give you further information with regard to its use.

Before using, please first read these instructions carefully.

What You Should Know About The SALIMA Bodywave® Sound Recliner/Mat

The SALIMA Bodywave® Sound Recliner/Mat combines the traditional sound bowl therapy with the latest creative, technological options.

Sound bowl therapy originated in the Himalaya region where it is traditionally practised in the field of healing and meditation. Sound bowls of various sizes are placed on the body and lightly tapped. The resulting vibrations are channelled through the body and in this way the oscillations penetrate deep into the tissue.

We have further developed the various options provided by this principle in the SALIMA Bodywave® Sound Recliner/Mat. Via a number of body sound generators integrated in the SALIMA Bodywave® Sound Mat there arises, in the low frequency range of 24-60Hz., an intensive cocoon like sensation. The exceptionally deep levels of penetration which result in this process ensure high rates of relaxation and harmonisation.

Six differing programs, each of 22 minutes duration are provided for selection. In addition to whole body application, individual body areas can be specifically addressed (the programs „Tension relief“, „Lower Back And Organs“ and „Breast And Neck“).

Apart from the relaxation element, further programs are available for revitalisation (the Programs „Vitalising“ and „Happy Day“ and which are designed to enhance energy levels to provide the basis for an active start in the new day.

General

Purpose Of Use

The SALIMA Bodywave® Sound Recliner/Mat with its pertinent operating parts has been designed to provide relaxation and regeneration for the user.

Safety Instructions

Keep water and other liquids well away from the operating part and electrical connections in order to avoid electric short circuiting.

Do not mishandle the SALIMA Bodywave® Sound Recliner/Mat and its operating part and do not expose to extreme temperatures nor high levels of humidity.

Areas Of Application (Indications)

- » Physical or psychological tranquility needs.
- » Physical or psychological tension or exhaustion.
- » Muscular tension.
- » After lengthy periods of sitting or standing.
- » Moodiness.
- » Stressful situations.
- » To underpin therapy, for example within the framework of pain relief treatment, physiotherapy or with infusions (on account of positive experience in these fields) or
- » simply for the purpose of relaxation and re-invigoration.

Contra-Indications

In respect of the following contra-indications, the SALIMA Bodywave®-Sound Recliner/Mat should only be used after consultation with your doctor.

- » Epilepsy and associated disorders (fits).
- » Wearers of electronic implants.
- » Serious mental illnesses.

Side Effects

There are no known undesirable side effects.

Description Of The SALIMA Bodywave®-Sound Recliner/Mat

1. Control Unit
2. Toggle Switch ON/OFF
3. Power Connection Socket
4. 15 pole Connection Socket For The SALIMA Bodywave® Sound Recliner/Mat
5. Rotary Volume Control Knob
6. LED Control Lamp Green
7. 3.5 mm Jack Socket For External Audio Sources

Black Variation

Beige Variation

Beige Sound Mat

Start-Up

Control Unit ①

The control unit is connected to the SALIMA Bodywave® Sound Recliner/ Sound Mat.


- » Insert the plug on the cable for the SALIMA Bodywave® Sound Recliner/Mat into the 15 pole socket ④ on the reverse side of the control unit and fasten with the screws. The plug is easily inserted into the socket in the one direction stipulated.

Power Connection

- » Insert the mains power plug into the power socket ③ on the control unit.
- » The mains cable is to be connected with the mains power supply socket (230 V/50Hz.).
- » Switch on the control unit with the ON/OFF switch on the reverse side of the casing. ②
- » green control lamp ⑥ right, down, on the front side of the control unit, will illuminate.

Master Display

- » Please wait until the master display with program selection appears.

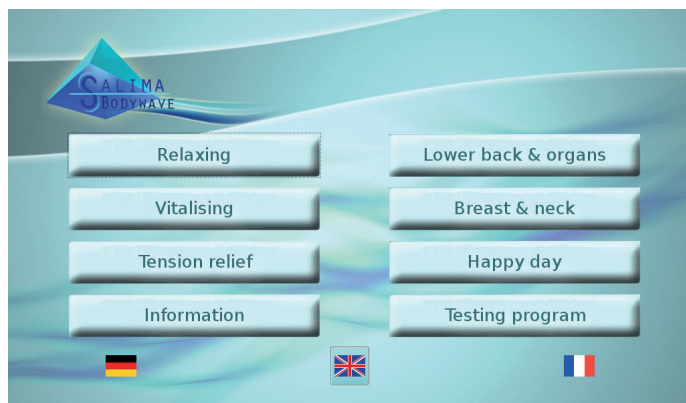Program Selection

- » Six programs for relaxation and revitalisation are available for selection.
- » The information program shows function and program sequence.
- » The test program serves to check the function of the body sound generator.
- » Depress the button for the selected program once to enable the program to start.
- » The system is now operational and the program is loaded. Please wait until the start screen with the program selection appears.

Recommendation: Run the test program on a regular basis in order to check the function of the body sound generator which is integrated in the sound recliner. The body sound generator highlighted in the sequence must vibrate simultaneously with the display.

If not all four body sound generators are perceptible please get in touch with your customer service facility.

Buttons and Controls

The sliding controller which sets the intensity of the vibration is located on the right hand side of the screen. By touching lightly with the finger tips the slide can be drawn upwards and downwards, (see illustration below).

- » By touching the „Pause“ button ⑧ the program goes on hold; re-newed touching re-starts the program.
- » By touching the „Stop“ button ⑨ the program which is running is halted and the program selection window re-appears.
- » With the sliding controller ⑩ on the right hand side of the display the vibration intensity can be set in accordance with the user's requirements.
- » To listen to one's own music on the SALIMA Bodywave® Sound Recliner/Mat, a 3.5 mm cable can be connected, via the jack socket, to an external audio source, for example an MP3 player, mobile telephone. Such cable is not included in the shipment.
- » After 22 minutes of running time, the starting window re-appears automatically.
- » Depress the toggle switch to switch off the control unit ② on the reverse side of the unit.
- » The volume controller ⑤ left, down, on the control unit is used to regulate the music volume level according to personal requirements.

Application Duration

The duration of the application program is 22 minutes and this should not be shortened.

Recommendation: When the program is completed, get up again slowly.

Programs Overview***Program 1 - Relaxation***

Constant pulsed vibrations from the buttocks to the apex, relaxing and soothing.

Program 2 - Revitalisation

Vibration on all parts of the body in rapid alternating sequence, revitalising and energising.

Program 3 - Relief For The Legs

For the whole back area, particularly for the legs, relaxation after long periods of sitting or standing. Can stimulate the lymphatic circulation, harmonising.

Program 4 - Lower Abdomen, Lumbar Vertebral Column

Vibration wave permeates the whole of the body, particularly in the region of the lower abdomen, cramp easing, blockades can be unblocked specially in the area of the lumbar vertebral column and the legs.

Program 5 - Neck/Breast

Vibration wave migrates through the whole of the body, particularly in the region of the upper body. Relaxes the cervical vertebral column and the thoracic spine. Mood compensating, relaxation promoting right down into the legs.

Program 6 - Happy Day

Rapid interchange of tangible vibration waves. Relaxing, harmonising, stimulating and invigorating.

General Information

Care And Cleaning

No special care and cleaning agents are required for the SALIMA Bodywave® Sound Recliner/Mat and the control unit. We recommend cleaning with, for example, a micro fibre cloth.

Care should be taken to ensure that no moisture penetrates the interior of the control unit. Should this nevertheless occur the unit must be subjected to a safety technical check before being re-used.

Delivery Content

1 SALIMA Bodywave® Sound Recliner complete system	REF-Nr. 200179
1 Mains connecting cable	REF-Nr. 200204
1 Neck Pillow	REF-Nr. 200202
1 Foot Rest Transparent	REF-Nr. 200200
1 Operating Instructions, multi language	REF-Nr. 101763
1 Abbreviated Operating Instructions, German	REF-Nr. 101741
Design and colour optional	
1 SALIMA Bodywave® Sound Recliner Beige	REF-Nr. 200179
1 Beige1 SALIMA Bodywave® Sound Recliner black	REF-Nr. 200179-black
1 SALIMA Bodywave® Sound Mat beige	REF-Nr. 200183-beige
1 SALIMA Bodywave® Sound Mat black	REF-Nr. 200183 black
1 SALIMA Bodywave® Sound Recliner Arm Rest Left	REF-Nr. 200210-L
1 SALIMA Bodywave® Sound Recliner Arm Rest Right	REF-Nr. 200210-R

Resettings, Modifications and Repairs

The manufacturer is only then responsible for the safety and performance of the device providing any resettings, modifications and repairs are carried out by personnel authorised by the manufacturer and that the product has been operated in accordance with the operating instructions.

Warranty

The normal warranty rights in accordance with BGB apply.

Guarantee

The manufacturer guarantees the SALIMA Bodywave® Sound Recliner/Mat for a period of twelve months from the date of acquisition by the end customer.

The guarantee is not applicable for wear and tear items such as the connecting cable, the recliner itself in case of damage which results from incorrect operation of which the customer was aware and damage resulting from the customer's own fault.

Technical Data

Nominal voltage	100 - 240 VAC
Nominal current	1.5 A
Operating frequency	50 - 60 Hz
Display	10" TFT Touch Display
Resolution	1024 x 600 Pixel
Number of body sound generators	4 pieces
Frequency range of the body sound generator	24 - 60 Hz
Mains fuse	T 2.5 A 250 V

Explanation Of The Illustrations On The Rating Plate

Attention: Read accompanying documents, especially user manual!

Follow instructions for use

Catalogue number

Serial number

Environment

Do not dispose of the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this you will help to protect the environment.

By labelling with CE mark, the manufacturer asserts that the product fulfills all relevant requirements of the respective EC Directives. A conformation assessment process has been successfully completed.

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Pierenkemper GmbH